



Scandi
orthopedic

Date

 | |

Measurement form DogLeggs Standard length adjustable

Veterinarian / therapist

Name

E-mail

Tel

Clinic

Name

Address

VAT-number

Patient

Name

Breed

☐ Dog

☐ Female

☐ Cat

☐ Male

Age

Weight

Diagnosis

Other medical issues

Measuring

Step 1 - Size calculation Standard length Adjustable Dogleggs
Measure as instructed on page two.


Step 2 - Size guide Stabilization system for Shoulder joint
Order Stabilization system for Shoulder joint based on the size for Standard length Adjustable Dogleggs, see page three.

Step 1

Measure with animal standing


Measured in ☐ cm ☐ inches

NOTE! For the breed english bulldog a custom order is needed.

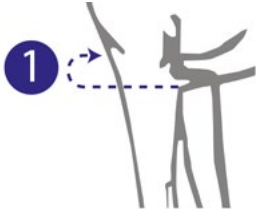


1. Measure the circumference of each leg at the point of the elbow. Use the largest of the two measurements. Round off to the closest integer.

Result #1

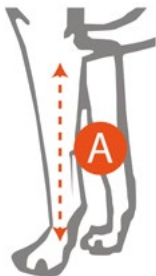

A. Measure the distance between the elbow and the top of the paw.
If this distance is:


2

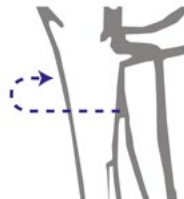

- **More** than 15 cm / 6 inches, then measure the circumference of each leg 10 cm / 4 inches below the point of the elbow. Round off to the closest integer.


Result #2

2

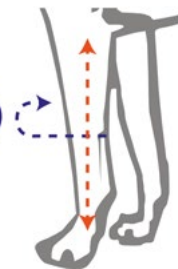

- **Less** than 15 cm / 6 inches, then measure the circumference of each leg at half the distance (A). Round off to the closest integer.

Result #2

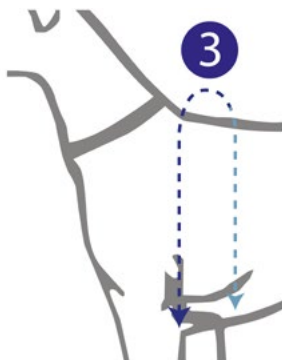

3. Measure from the point of elbow on one side, over the back, to point of elbow on the opposite side. Divide this number by two, then add two. ($x/2 + 2$). If the sum is an odd number then round up to the closest even number (ex 4,6,8 etc).

Result #3

Apply the results like this: Result #1 / Result #2 / Result #3

Size:

 / /

Stabilization system for shoulder joint

The Stabilization system can only be used together with Standard length Adjustable Dogleggs.

If size for Std length adjustable begins with:	Then the right size for Stabilization system is:
1 - 4	18.05 ☐
5	18.1 ☐
6	20.2 ☐
7	28.2 ☐
8	32.2 ☐
9	34.2 ☐
10	36.2 ☐

Scandi Orthopedic AB

Christoffer Diehl

Certified Prosthetist/Orthotist (CPO) and CEO

Address

Scandi Orthopedic
Ringögatan 5
417 07 Göteborg

Telephone: 031-21 41 28

E-mail: info@scandiop.se

www.scandiop.se

© 2020 Scandi Orthopedic AB

